

Procedure Guidelines for Transoral Removal of Supraglottic Laryngeal Tumors using Microelectrodes

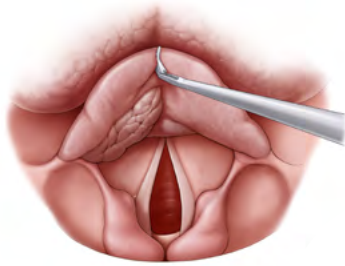


Fig. 1: Supraglottic carcinomas

Indications and contraindications

The indications and contraindications for resections of supraglottic laryngeal tumors using ARROWtip™ monopolar microdissection electrodes correspond to those for the surgical removal of benign and malignant supraglottic tumors. However, it is a prerequisite that the patient can be properly prepared for resection of the tumor with the micro-laryngoscope or distending laryngoscope. Tumors on the exposed edge of the epiglottis, central suprahoid epiglottic carcinomas, carcinomas at the edge of the false vocal cords, and carcinomas in the area of the aryepiglottic fold are especially well suited for this (Fig. 1).

Patient preparation

General anesthesia is used. A local anesthetic (e.g. Ultracaine 1% Suprarenine) may be injected with a vasoconstrictor, depending on the surgeon's preferences. Since monopolar electrodes are used, a neutral electrode is attached to the patient (e.g. right upper arm).

Intervention

The surgical procedure begins with an operating laryngoscope (e.g. Kleinsasser micro-laryngoscope). The hand instruments still required are micro-forceps and a non-stick suction tube (REF 71 50 17) for hemostasis, which are also used in conventional laryngeal microsurgery.

The benefit of ARROWtip™ monopolar microdissection electrodes is that resection outside the line of sight is also possible thanks to the available angled models.

Smaller, defined epiglottic carcinomas (suprahoid) can be resected using an ARROWtip™ (straight or angled; REF 36 44 71 or REF 36 44 75) along the tumor with the required distance.

The depth extension of epiglottic carcinomas (infrahyoid) must be considered prior to the procedure. This is often difficult in the area of the petiole. The preepiglottic space may be infiltrated with such carcinomas.

In this case the preepiglottic adipose body first has to be exposed after severing the vascular bundle in the pharyngoepiglottic fold with Calvian® bipolar forceps (REF 70 09 46). This is performed by means of medial transection of the glossoepiglottic fold with a straight ARROWtip™.

In the next step, the epiglottis is transected with a straight ARROWtip™ monopolar microdissection electrode (Fig. 2). The epiglottis is pulled in the dorsal direction (towards the pharynx) using forceps during this step. This also allows the lingual surface to be seen.

Following the caudal sagittal incision of the epiglottis, a horizontal lateral incision is made in the next step (Fig. 3) with a 90° angled ARROWtip™ monopolar microdissection electrode (REF 36 44 73).

A good view of the incision is important here in order to avoid injuries to other structures (cartilage, nerves, larger vessels). If areas of the thyroid or arytenoid cartilage are also infiltrated, they can be included in the resection. A temporary precautionary tracheotomy and feeding by stomach tube may be required. This needs to be decided on a case-by-case basis.

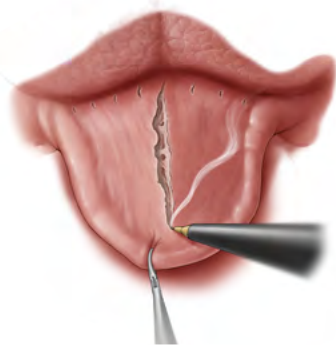


Fig. 2: Medial transection of the glossoepiglottic fold with a straight ARROWtip™ monopolar microdissection electrode (REF: 36 44 71).


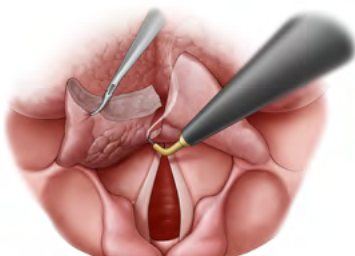


Fig. 3: Horizontal lateral incision with a 90° angled ARROWtip™ monopolar microdissection electrode (REF: 36 44 73).

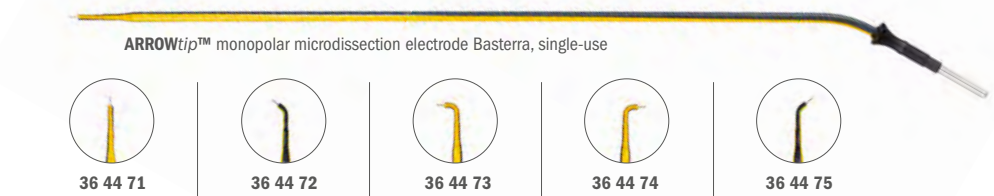

Postoperative treatment

The patient is discharged 1 to 5 days postoperatively depending on the extent of resection. A temporary precautionary tracheotomy should be performed for larger resections.


Partial or complete voice rest in the first few days is recommended. Logopedic exercises may already be indicated at this stage as well. Antitussives may be used to suppress a nervous cough and corticosteroids if there is a risk of edema. Postoperative treatment is similar to that of patients who undergo the CO2 laser procedure.

Disclaimer: These procedure guidelines have been carefully researched and compiled with the help of specialist physicians. They are not meant to serve as a detailed treatment guide. They do not replace the user instructions for the medical devices used. Sutter accepts no liability for the treatment results beyond legal regulations.

Settings* for CURIS® 4 MHz radiofrequency generator (REF: 36 01 00-01)



ARROWtip™: Monopolar CUT 2
Power adjustment: 5 to 25 watts

Calvian®: Bipolar MACRO
Power adjustment: 20 watts




For further accessories see back page.

Valid for the CURIS® with the orange label.

ARROWtip™: Monopolar CUT 2
Power adjustment: 25 to 46 watts


Calvian®: Bipolar MACRO
Power adjustment: 20 watts



For further accessories see back page.

* Always start with the lowest settings to achieve the desired effects. If necessary, increase the settings step-by-step until the desired effect is achieved. This may even be 50 watts or higher. The settings may differ from patient to patient, from tissue to tissue, and have to be adjusted accordingly.
Please consider that this information is not meant to serve as a detailed treatment guide.

Recommended products for this treatment

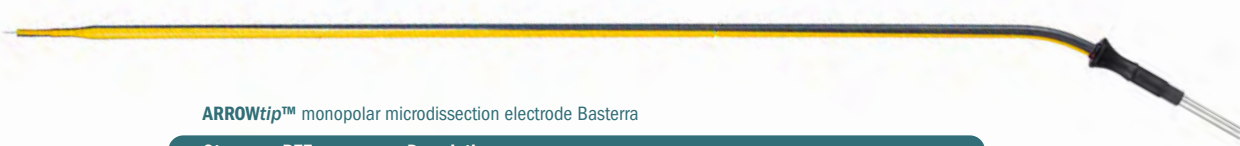

ARROWtip™ monopolar microdissection electrode Basterra

Qty.	REF	Description
2	36 44 71	ARROWtip™ monopolar microdissection electrode Basterra, single-use working length: 212.0 mm, total length: 232.0, straight tips
2	36 44 72	ARROWtip™ monopolar microdissection electrode Basterra, single-use working length: 210.0 mm, total length: 229.0, 45° angled downwards, Ø 2.4 mm
2	36 44 73	ARROWtip™ monopolar microdissection electrode Basterra, single-use working length: 207.0 mm, total length: 227.0, 90° angled downwards, Ø 2.4 mm
2	36 44 74	ARROWtip™ monopolar microdissection electrode Basterra, single-use working length: 207.0 mm, total length: 230.0, 90° angled upwards, Ø 2.4 mm
2	36 44 75	ARROWtip™ monopolar microdissection electrode Basterra, single-use working length: 210.0 mm, total length: 233.0, 45° angled upwards, Ø 2.4 mm



Calvian® bipolar clamp

134° C
autoklavierbar




Qty.	REF	Description
1	70 09 46	Calvian® bipolar clamp, working length: 23.0 cm, 45° angled tips, with suction



CURIS® 4 MHz radiofrequency generator
Basic Equipment

Qty.	REF	Description
1	36 01 00-01	CURIS® 4 MHz radiofrequency generator (incl. mains cord, user's manual and test protocol)
1	36 01 10	Foot switch two pedals for CURIS® (cut & coag), 4 m cable
1	37 01 54 L	Bipolar cable for CURIS®, length: 3 m
1	36 07 04	Monopolar handpiece (pencil) cut & coag, shaft 2.4 mm, cable 3 m
1	36 02 38	Cable for single-use patient plates, length: 3 m



available patient plates:

1 (x100)	29 00-5	Single-use patient plate, split, for adults and children, PU 20 x 5 pcs.
1 (x50)	95 80 04	Single-use patient plate, split, for adults, PU 10 x 5 pcs.
1 (x50)	95 80 05	Single-use patient plate, split, for children, PU 10 x 5 pcs.
1	36 02 26	Re-usable rubber patient plate

Product availability is subject to regulatory approval in individual markets. Products may therefore not be available in all markets. Lengths for orientation purposes; may vary slightly.



PRECISION ELECTROSURGERY
Made in Germany

SUTTER MEDIZINTECHNIK GMBH

TULLASTRASSE 87 · 79108 FREIBURG/GERMANY
TEL. +49(0)761-51551-0 · FAX +49(0)761-51551-30
WWW.SUTTER-MED.COM · INFO@SUTTER-MED.DE